

# AVChoice™

DEDICATED VETERINARY IMAGING SYSTEM  
High Frequency, Anatomically Programmed



AVChoice systems are equipped with an integrated tube stand and an anatomically programmed, high frequency generator providing a cost effective and time saving solution for the veterinarian who seeks maximum capability in minimum space.



## Key Features

- Integrated tube stand with variable SID
- Angulating tube arm
- Control console on swing out arm mounts to tube stand
- Available with stationary, 2-way float or 4-way float table top

Exposure footswitch

# AVChoice™

DEDICATED VETERINARY IMAGING SYSTEM



The operator control panel with APR (Anatomical Programmed Radiography) technique setup conveniently mounts on a swing-out arm off either side of the tube stand. This gives the operator the ability to bring the control panel forward or position it to the side for convenient technique setup.

## Standard Features

- Welded table construction with 400lb load capacity and urine trap
- Integrated tube stand with electric locks and variable SID travels full length of table
- Angulating tube arm, angulation dial and operator handle
- Electric table locks (on 2-way float and 4-way float table tops)
- Foot activated exposure switch
- 17in x 17in grid cabinet with tube stand interlock
- Grid: 10:1 ratio, 103 line, 34-44in focal range
- Heavy-duty manual cassette tray
- Certified manual collimator with built-in spare bulb and retractable tape measure
- X-ray tube: 140,000 heat unit, 1.0-2.0mm focal spots with 90°horn angle
- One pair 20ft high voltage cables, terminated both ends

## System Options

- Available with stationary, 2-way float or 4-way float table top
- DelWorks flat panel DR upgrade package


## Generator

- ANTHEM 500mA, 32kW, high frequency, three-phase with anatomical programming


## Generator Options

- Pedestal stand or wall mount for console
- Electronics mount in base of table
- Single phase and 40kW generator options also available


AVChoice all models

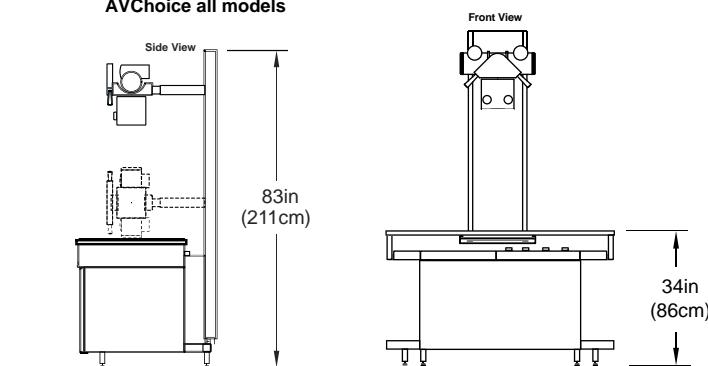

AVChoice stationary and two-way float

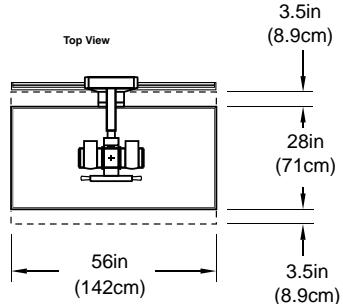

AVChoice four-way float

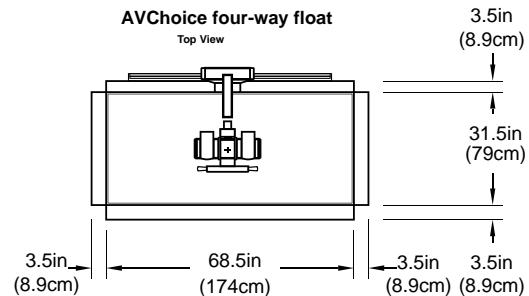

## DEL MEDICAL

Toll Free: (800) 800-6006  
Fax: (847) 288-7011  
Email: [info@delmedical.com](mailto:info@delmedical.com)  
[www.delmedical.com](http://www.delmedical.com)  
241 Covington Drive, Bloomingdale, IL 60108

## Universal

Universal branded products are manufactured and distributed by UMG/Del Medical. Del Medical is an ISO 13485 & ISO 9001 Certified facility. Specifications are subject to change without notice.

## UMG

Toll Free: (800) 261-9808  
Fax: (914) 835-6111  
Email: [info@umgxray.com](mailto:info@umgxray.com)  
[www.umgxray.com](http://www.umgxray.com)  
28 Calvert Street, Harrison, NY 10528